



Specification

Model	VS-10 S, M, H
Dimensions	191(H) x 92 (W) x 112 (D) mm
Weight	About 0.25kg (including the battery)
Power Adapter	Input voltage: AC 100-240V 50/60Hz Input power: 25VA. DC output voltage: 5V 2A
Battery	Built-in battery: 3.6V 3400mAh Continuous operation duration of the lithium battery: 4 hours Lithium battery charge time: not longer than 4 hours
Rotation Angle of Display Screen	Maximum vertical rotation angle: 140°±10° Maximum horizontal rotation angle: 270°±10°
Screen	3.5 inches touchscreen with resolution of 640x960 pixels
Water Proof Level	IP66
Photo Record	Photo can be recorded in JPEG format, up to 10,000 photos can be recorded
Video Record	Video can be recorded in MPEG-4 format, up to 2 hours video can be recorded
Media Storage	Not less than 5GB
Data Port	USB Type C, HDMI ^① , and Wi-Fi ^② are available
Anti-Fog	Immediate anti-fog technology is supported with no need for heating of the anti-fog coating on the laryngoscope blade
Blades	Single use blade
Frequency for WiFi ^②	2.4GHz

Note: ① For VS-10H & VS-10M, ② For VS-10H only

Smaller Size, Bigger Vision
VS-10 Series Video Laryngoscope*



MEDCAPTAIN MEDICAL TECHNOLOGY CO., LTD.



Add.: 12th Floor, Baiwang Research Building, 5158 West Shahe Road, Xili, Nanshan, 518055 Shenzhen, Guangdong, PEOPLE'S REPUBLIC OF CHINA
Tel: +86 755 26953369 E-mail: info@medcaptain.com
Website: www.medcaptain.com
MEDCAPTAIN is a trade mark of Medcaptain Medical Technology Co., Ltd.
Specifications are subject to change without prior notice.
©2023 Medcaptain Medical Technology Co., Ltd. All rights reserved.
EN-VS-10 Series-4P-Version 1.6

Distributor:



VS-10 Series Video Laryngoscope*



reddot winner 2021



HDMI: Bigger screen display

Vision with Safety



Smart Power Management
Remaining Time Display



IP66, Waterproof Design



Instant anti-fog technology
No heating time required

Vision with Ease



3.5"HD touchscreen: a step
closer to reality



Horizontal rotation angle: 0°~270°



Vertical rotation angle: 0°~140°



Wifi: Data transmission and projection
to bigger screen

Vision For Everyone

All 5 blades share a general handle
Disposable blades offer a better infection control



Size **M1** for infant



Size **M2** for child




Size **M3** for adult



Size **M4** for obesity/big adult



Size **M3D** for difficult intubation

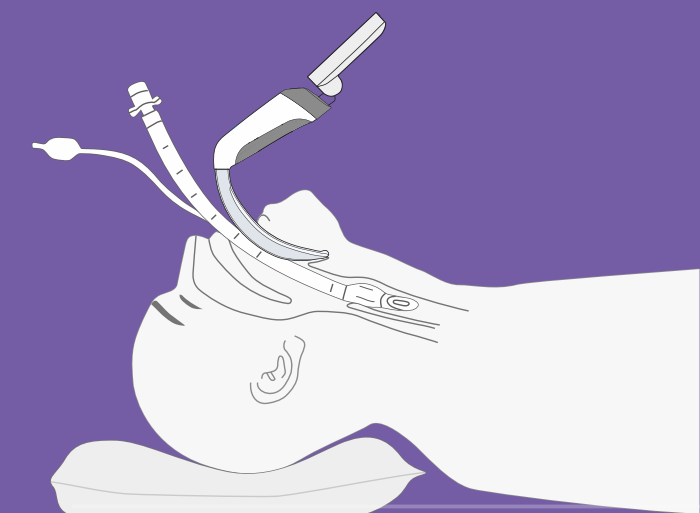

*Not available in Germany